Preparing the Scale for Use


- 1. Remove the battery cover on the back of the scale and insert a 9V battery (Not included) . Replace cover.
- 2. If you use the scale to weigh an infant who cannot yet stand on his own, attach the tray to the scale. Then turn the locking nut counter-clockwise to secure the tray in place.

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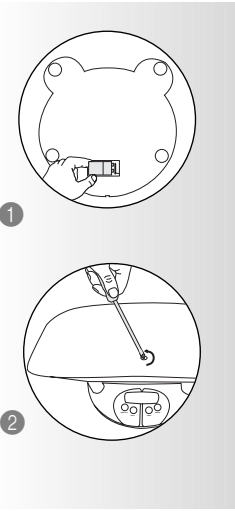
Operating the Scale

Note: Place the scale on a hard, flat surface. Do not place the scale on a rug or soft surface.

- 1. Press “ON/OFF”. The display shows “8888”.
- 2. Wait until the display shows “0:0” or “0.00”.

Note: If the display shows “0:0”, the scale is set to pounds/ounces. If the display shows “0.00”, the scale is set to kilograms/grams.

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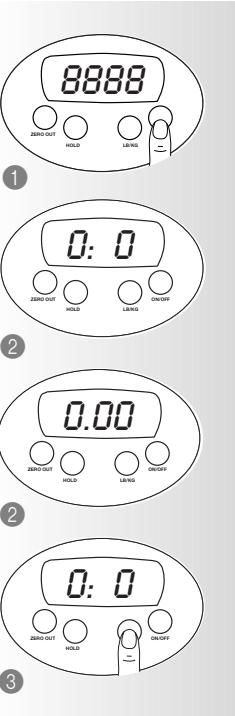
Using the Scale With the Tray

Caution: Before use, be sure that the tray is securely locked in place.

- 1. Press the “ON/OFF” button and follow the instructions in the “Operating the Scale” section.
 - 2. Lay the baby on the tray and read the display. If the baby is moving, press the “HOLD” button to lock in the correct weight. (See “Using the ‘Hold’ Feature.”)
- Note: The scale will automatically turn off after 2 minutes.

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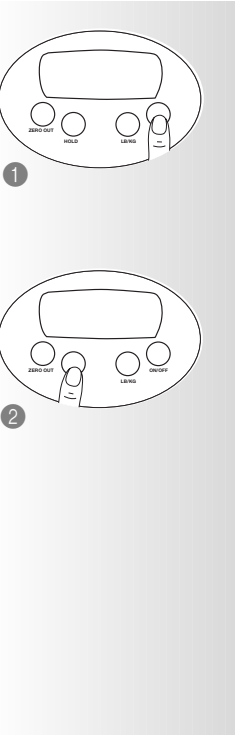
Using the Scale Without the Tray

- 1. Remove the tray by turning the locking nut clockwise.
- 2. Insert the 4 rubber plugs and insert the supplied plastic cover tab on the locking tray hole. Follow the instructions in the “Operating the Scale” section.
- 3. Hold your baby upright and place the baby's feet so that the weight is evenly distributed. The LCD display will display his/her weight.
- 4. If necessary, press the “HOLD” button to lock your baby's weight.

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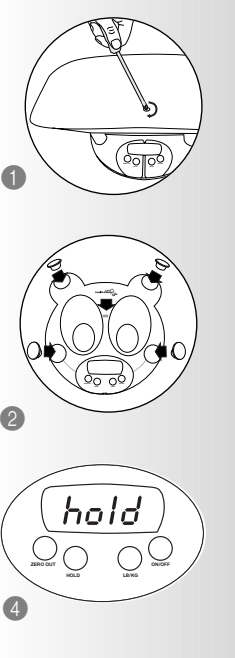
Using the “Zero Out” Feature

If you are weighing the baby with a blanket or diaper:

- 1. Place the blanket or diaper on the tray.
- 2. Press the “ZERO OUT” button until the display reads “0:0” or “0.00”.
- 3. You can now accurately weigh your baby.

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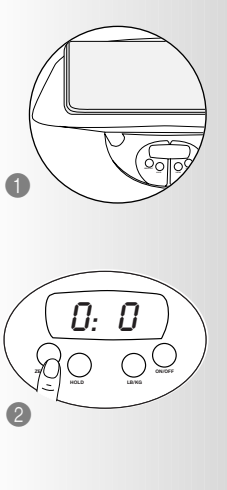
Using the “Hold” Feature

If your baby is moving:

- 1. Press the “HOLD” button.
- 2. The display will first read “hold,” then the baby's average weight will display. (The “▲” symbol will flash on the display.)

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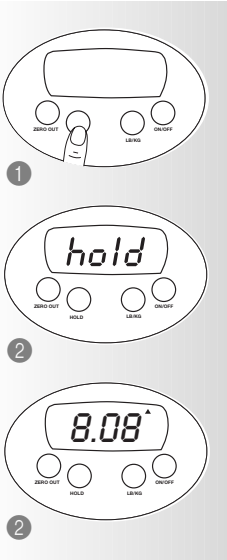
Taking Care of Your Scale

- 1. Periodically wipe your scale with a damp cloth to remove dust.
- 2. Be careful handling the scale. It is a delicate instrument.
- 3. Do not use abrasive cleansers. They will damage the scale's surface.

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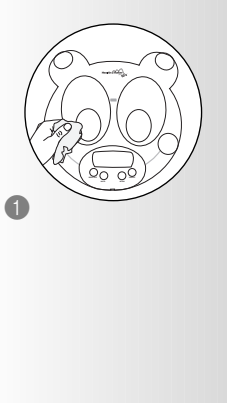
Symbols

- “Lo” appears when the battery needs to be replaced.
- “Err” appears when there is too much weight on the scale.
- If “---” appears, turn off the scale and turn it back on. Wait until “0:0” or “0.00” appears to start weighing.

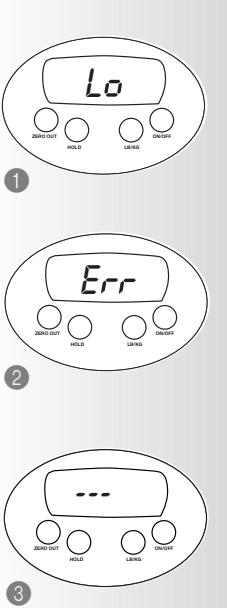Caution:

Weight exceeding the capacity (60lb/27kg) may damage your scale. This scale is designed for household use only. This scale is not waterproof; avoid contact with excessive moisture.

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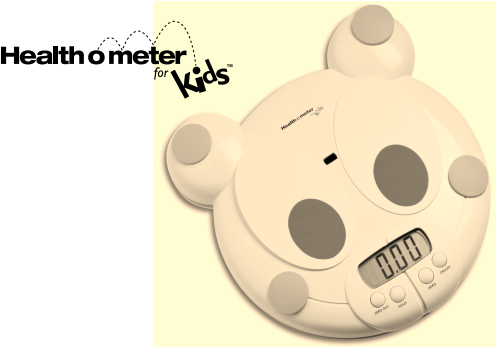
READ ALL INSTRUCTIONS BEFORE USING THIS PRODUCT.

SAFETY INFORMATION

- The Health o meter for Kids™ scale is intended for home use with adult supervision. Never leave a child unattended on the scale.
- Do not place the unit on slippery surfaces.
- Ensure the tray is properly locked before placing the baby on the scale.
- Do not immerse into water or pour liquids onto the scale.
- Do not use while bathing.

Note: It is normal for weight to vary during the day and from one day to another. Slight changes in weight are generally normal. If you are concerned with your child's weight fluctuation, please contact your physician.

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Thank you for Buying the Health o meter for Kids™ scale!

Congratulations on your purchase of the Health o meter for Kids™ HDC100 Baby/Toddler Scale. This special scale is designed to accurately and reliably weigh your baby in both pounds and kilograms, up to a maximum of 60 lbs (27kg.)

- You can wrap your baby in a soft blanket and still get an accurate reading. The scale lets you easily deduct any added weight with the “Zero Out” feature.
- Even if your baby moves, the scale's “Hold” feature displays a steady, average weight.
- You can easily remove the holding tray when your baby gets older and begins to stand.

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© 2003 Sunbeam Products, Inc. All rights reserved. Health o meter for Kids™ is a trademark of Sunbeam Products, Inc.