
# Health o meter® BODY FAT MONITORING SCALE

User's Manual  
Model BFM580
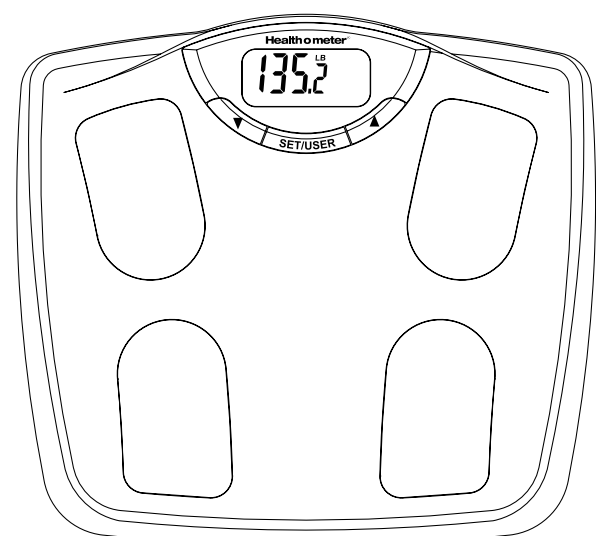
### PLEASE READ CAREFULLY AND SAVE


For information on service or warranty, contact Sunbeam Products, Inc. at 1-800-672-5625 or go to [www.healthometer.com](http://www.healthometer.com).

To assist us in serving you, please provide model number and date of purchase.

## Monitoring your Weight and Body Fat Composition

\* Make sure you have programmed your personal data.

\* Place the scale on a flat, level surface. A hard surface is best to achieve maximum accuracy. Placing the scale on a rug or carpet may affect the scale's weighing accuracy.

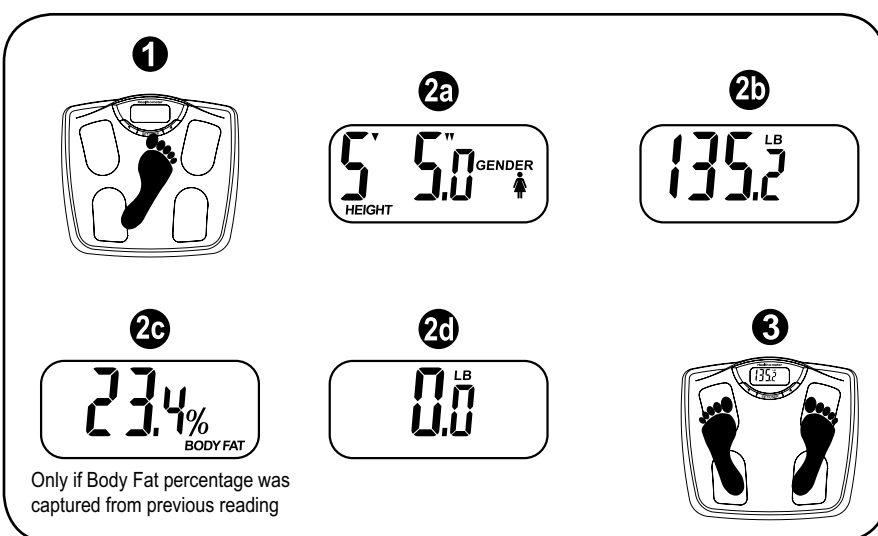1. Press the **SET/USER** button with your toe. Your personal profile will appear on the screen.

2. Wait until "0.0" appears on the display.

3. Stand on the scale (bare feet only!) so the balls of your feet are on the upper silver pads, heels are on the lower silver pads.

**NOTE: For accurate body fat estimates, always measure with bare feet.**

4. During the weight measurement, the display will flash "--" and your weight and body fat percentage will be shown. Stand still for a few seconds.

**NOTE: It is normal for your weight to vary during a day and from one day to another, slight changes in weight are generally normal.**

4

## WELCOME

Congratulations on your purchase of this Health o meter® wellness monitoring product. It's so much more than just a scale. It helps monitor the weight and body fat composition for one user. It's easy to program and use. Once your data is stored, your scale is ready to use. It's that simple!

## ABOUT BODY FAT

The Health o meter® Body Fat Monitoring Scale can help you in achieving your physical best by providing two critical measurement estimates: **BODY FAT COMPOSITION** and **BODY WEIGHT**.

Traditional bathroom scales, indicate your overall weight—both lean mass and body fat together. The Health o meter® Body Fat Monitoring Scale lets you track your weight and body fat separately.

### IMPORTANT!

Please read all instructions to familiarize yourself with features and operations before programming the scale.

### CAUTION!

- The use of this scale is not recommended for persons with an electric implant (heart pacemaker).
- The scale will not calculate body fat unless you have bare feet.
- For an accurate body fat estimate, do not bend your knees, and keep your legs/thighs apart and straight when standing on the scale.
- The body fat function is not recommended for use by pregnant women, persons with fever, swollen legs or other edemas, as well as a person who is over hydrated or dehydrated.
- Body fat percentages may be slightly higher or lower for children under 16, persons with diabetes and other health conditions.
- Women naturally store approximately 5% more body fat than men.

1

## Weight-only Usage

It is not necessary to have the user's personal data programmed into the scale to obtain an accurate weight-only reading.

\* Place the scale on a flat, level surface. A hard surface is best for maximum accuracy. Placing the scale on a rug or carpet may affect the scale's weighing accuracy.

- Tap on the scale platform with the ball of your foot to activate the scale.
- Wait until "0.0" is displayed.
- Step on the scale and stand still. The display will flash "--" and your weight will be shown.
- The scale will turn off automatically after a few seconds.

**NOTE: It is recommended to weigh yourself at the same time of the day. It is normal for your weight to vary during a day and from one day to another, slight changes in weight are generally normal.**

### CAUTION!

If "E" appears on the display, an error has occurred. Step off the scale and carefully repeat the steps again. The "E" will appear if:

- You step on the scale before "0.0" displays.
- The weight is greater than the maximum capacity (330LB / 150kg).  
NEVER apply a weight beyond the maximum capacity.
- The scale cannot measure body fat.

## Caring for your Scale

Periodically wipe the scale surface down with a DRY cloth to remove dust.


5

## OPERATING INSTRUCTIONS


### Installing the Batteries

- Locate the battery cover on the back of the scale.
- Remove the battery cover and insert the 4 AA batteries.
- If "Lo" appears in the display, replace the batteries.

**Note: Do not use rechargeable batteries.**
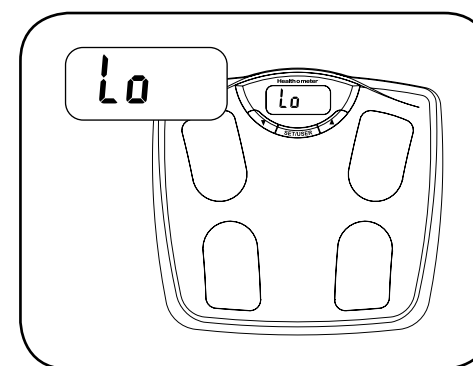
### Keys and Icons

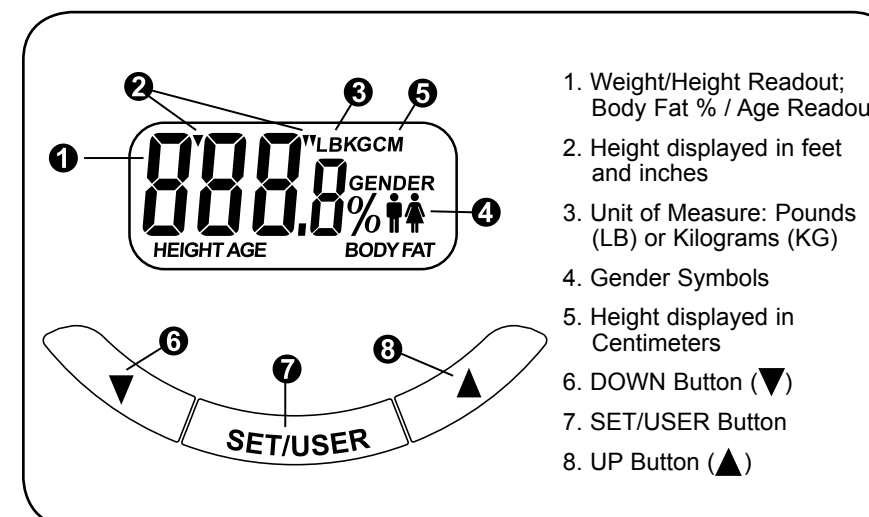
- Weight/Height Readout; Body Fat % / Age Readout
- Height displayed in feet and inches
- Unit of Measure: Pounds (LB) or Kilograms (KG)
- Gender Symbols
- Height displayed in Centimeters
- DOWN Button (▼)
- SET/USER Button
- UP Button (▲)



2

## Body Fat Composition Ranges

Men's Body Fat Composition Ranges

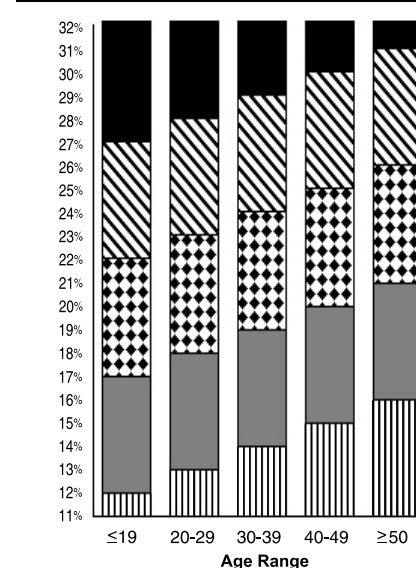Women's Body Fat Composition Ranges

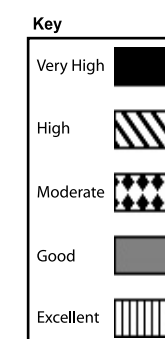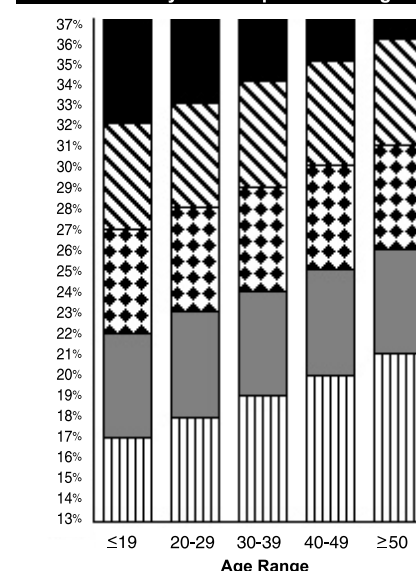Extracted from Principals & Labs for Fitness and Wellness, 1st Edition ©1999. To determine the percentage of body fat that is appropriate for your body, consult your physician.


6

## First Time Programming Instructions

On the back of the scale platform, use the selector switch to set the unit of measurement to "LB" (pounds) or "KG" (kilograms).

Turn the unit over and follow these steps to program your Body Fat Monitoring Scale:

- Press and hold the **SET/USER** Button until the **HEIGHT** icon starts flashing.
- Press the **▲** or the **▼** button to select the desired **FEET** (or **CENTIMETERS** in **KG** mode); then press the **SET/USER** button.
- Press the **▲** or the **▼** button to select the desired **INCHES**; then press the **SET/USER** button.
- Press the **▲** or the **▼** button until your correct **AGE** is shown on the display; then press the **SET/USER** button.
- Press the **▲** or the **▼** button to select your gender (**MALE** or **FEMALE**); then press the **SET/USER** button.
- Once you have entered your information, the scale will display your profile. The display will flash "--". Wait until the display shows "0.0". The scale is ready for use.

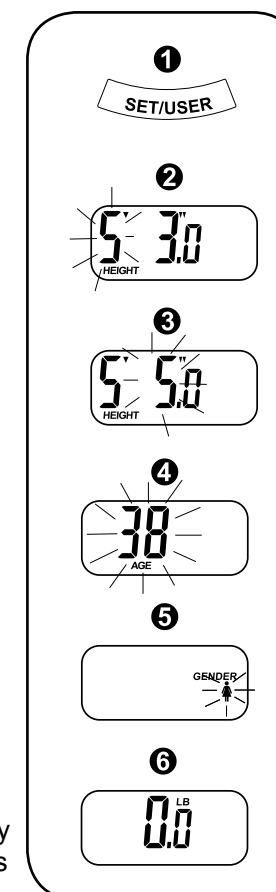**NOTE: If a button is not pushed within 20 seconds during programming, "0.0" will be displayed. Wait until the unit turns off and return to Step 1 to continue set up.**

**If the scale is operating in the Kilograms (KG) mode, the height will be set in centimeters (CM).**

3

## 5-YEAR LIMITED WARRANTY

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.**

Sunbeam Products, Inc. ("Sunbeam") warrants that for a period of five years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or in any way change the terms and conditions of this warranty. This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

### What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition.

Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

### How to Obtain Warranty Service

If you have any question regarding this warranty or would like to obtain warranty service, please call 1-800-672-5625, or visit us at [www.healthometer.com](http://www.healthometer.com) and a convenient service center address will be provided to you.

In the U.S.A., this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.

**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THE ADDRESSES SHOWN OR TO THE PLACE OF PURCHASE.**

©2005 Sunbeam Products, Inc. All Rights Reserved. Health o meter® is a registered trademark of Sunbeam Products, Inc., Boca Raton, FL 33431

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En E.E.U.U., esta garantía es ofrecida por Sunbeam Products, Inc. situada en Boca Raton, Florida 33431.