## FlexTemp™ Hot + Cold Joint Wrap

Model: 901-825-002

## IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY


**DANGER:** TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK AND FIRE, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS:

1. READ ALL INSTRUCTIONS.
2. DO NOT USE WHILE SLEEPING.
3. DO NOT USE ON AN INFANT OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, SLEEPING, HELPLESS OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, A PERSON WITH POOR BLOOD CIRCULATION OR INSENSITIVITY TO HEAT, OR A PERSON INCAPACITATED BY MEDICATION.
5. THIS APPLIANCE IS NOT TO BE USED BY CHILDREN. CLOSE SUPERVISION IS NECESSARY WHEN ANY APPLIANCE IS USED NEAR CHILDREN. CHILDREN SHOULD NOT PLAY WITH THE APPLIANCE.
6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.

7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING; CHECK SKIN UNDER PAD FREQUENTLY.
8. DO NOT USE IN AN OXYGEN ATMOSPHERE.
9. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
10. DO NOT SIT ON, AGAINST, OR CRUSH PAD – AVOID SHARP FOLDS. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
11. UNPLUG WHEN NOT IN USE. NEVER LEAVE APPLIANCE UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
12. NEVER PULL THIS PAD BY THE SUPPLY CORD. DO NOT USE THE CORD AS A HANDLE. DO NOT PINCH CORD.
13. DO NOT TUCK IN, TRAP, FOLD, CROSS, PINCH OR WRAP CONTROL OR CORD WITH HEATING PAD DURING USE.
14. USE CAUTION WITH PAD ENDS TO PREVENT HOOK FASTENERS FROM DAMAGING CLOTHING OR OTHER PROPERTY.
15. DO NOT PLACE ANYTHING OVER THE HEATING PAD WHILE IN USE OR PLUGGED IN, SUCH AS A TOWEL, BLANKET, CLOTHING, OR ANY OTHER INSULATING MATERIAL, AND NEVER ALLOW HEATING PAD TO WRAP AROUND ITSELF.
16. CAREFULLY EXAMINE PAD, CORD AND CONTROLLER BEFORE EACH USE. DISCARD THE PAD IF ANY SIGN OF DETERIORATION, WEAR OR DAMAGE. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CUSTOMER CONTACT INFORMATION.
17. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS. ALLOW TO COOL BEFORE STORING AND DO NOT CREASE BY PLACING ITEMS ON TOP OF IT. STORE IN THE ORIGINAL PACKAGING, IN A COOL, DARK, DRY PLACE.
18. IF PROVIDED WITH A DETACHABLE CONTROLLER, USE ONLY WITH THE PROVIDED CONTROLLER.
19. USE THIS PAD ONLY AT RATED VOLTAGE, AS MARKED ON THE PRODUCT. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.
20. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS. SKIN BURNS COULD RESULT.
21. THIS APPLIANCE IS NOT INTENDED FOR MEDICAL USE IN HOSPITALS.
22. SAVE THESE INSTRUCTIONS.

## HOUSEHOLD USE ONLY

### North American models with polarized plugs:

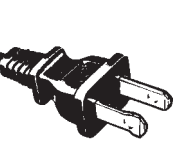**POLARIZED PLUG**  
This appliance has a polarized plug, (one blade is wider than the other). As a safety feature to reduce the risk of electrical shock, this plug is intended to fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature or modify the plug in any way. If the plug fits loosely into the AC outlet or if the AC outlet feels warm do not use that outlet.

## ADDITIONAL WARNINGS FOR SUNBEAM GEL PACK

1. READ AND FOLLOW ALL DIRECTIONS CAREFULLY TO AVOID POSSIBLE INJURY.
2. CONSULT A PHYSICIAN PRIOR TO USE.
3. PEOPLE WITH POOR CIRCULATION, DIABETES, NERVE DAMAGE, PARALYSIS OR SENSITIVE SKIN SHOULD NOT USE THIS PRODUCT.
4. ALTHOUGH NON-TOXIC, DISCARD GEL PACK IF PRODUCT LEAKS OR IS PUNCTURED.
5. ADULT SUPERVISION IS REQUIRED FOR USE ON INVALIDS OR CHILDREN.
6. DO NOT USE ON INFANTS OR CHILDREN UNDER THE AGE OF 2 YEARS. NOT FOR USE IN CRIBS OR PLAYPENS.
7. DO NOT USE THE PACK IN CONJUNCTION WITH LINIMENTS, SPORTS CREAMS, BALMS OR LOTIONS.
8. DO NOT USE ON BROKEN SKIN.
9. DO NOT FALL ASLEEP WITH THE GEL PACK.
10. DO NOT INGEST OR INHALE THE GEL. IF INGESTED, SEEK MEDICAL ATTENTION IMMEDIATELY.
11. DO NOT USE GEL PACK WITH MESH POUCH WHEN THE PAD IS CONNECTED TO ELECTRICAL SOURCE. GEL PACK IS ONLY TO BE HEATED BY MICROWAVE.

The FlexTemp™ Hot + Cold Joint Wrap is specifically designed for use on the elbow or knee joints for hands-free, targeted hot or cold therapy. Suggested treatment duration for hot or cold application to the joint is 20 minutes, up to three times a day, unless otherwise directed.  
Note: Always consult with your physician or health professional for recommended treatment of an injury or chronic muscle or joint pain.

### HOT VS. COLD THERAPY

| Heat Therapy helps to:                                                                                                                                                        | Cold Therapy helps to:                                                                                                                                     |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <ul style="list-style-type: none"><li>• Relax tight muscles</li><li>• Ease stiffness</li><li>• Facilitate stretching</li><li>• Relieve common muscle and joint pain</li></ul> | <ul style="list-style-type: none"><li>• Reduce inflammation</li><li>• Ease pain from acute injuries such as strains, sprains, bumps, and bruises</li></ul> |

### JOINT WRAP COMPONENTS

Figure 1  
901-825-002

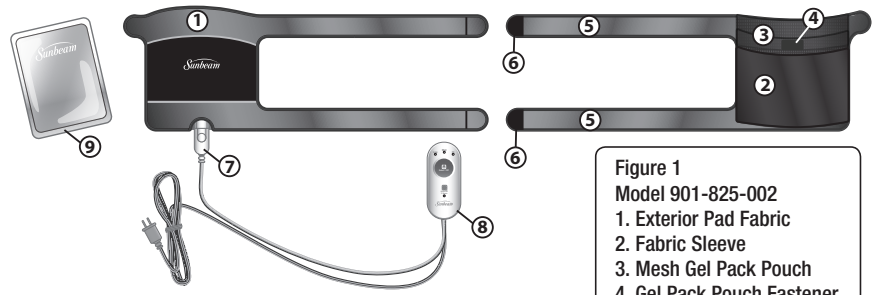
Figure 1  
Model 901-825-002  
1. Exterior Pad Fabric  
2. Fabric Sleeve  
3. Mesh Gel Pack Pouch  
4. Gel Pack Pouch Fastener  
5. Elastic Fastening Straps  
6. Strap Hook Fastener  
7. Controller Connector  
8. Controller  
9. Gel Pack


## USAGE INSTRUCTIONS

### HOW TO APPLY YOUR JOINT WRAP

Note: Some wear of the exterior pad fabric may be seen over time; this is normal and will not affect the performance of the heating pad. Use caution when fastening straps to prevent hook fasteners from touching clothing.


### HOW TO APPLY TO ELBOW – SEE FIGURE 2

1. Position wrap so Controller Connector (Figure 1, number 7) is hanging down. Insert hand into Fabric Sleeve (Figure 1, number 2) and pull the wrap up aligning the Sunbeam logo to the elbow. See Figure 2A.
2. The Elastic Fastening Straps (Figure 1, number 5) stretch for size adjustability and a comfortable fit. Use available hand to wrap each strap around the arm and secure in place by attaching the Strap Hook Fastener (Figure 1, number 6) directly to the Exterior Pad Fabric (Figure 1, number 1) or onto the strap itself. See Figure 2B.
3. To remove wrap, detach and unwrap Elastic Fastening Straps and pull arm out of Fabric Sleeve.
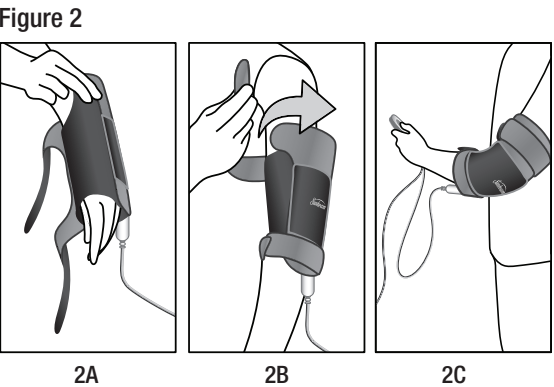
### HOW TO APPLY TO KNEE – SEE FIGURE 3

1. Position the wrap over the knee so the Sunbeam logo is aligned to the knee cap. The Fabric Sleeve (Figure 1, number 2) is not required for application to the knee.
2. Wrap each of the Elastic Fastening Straps (Figure 1, number 5) around the leg and secure in place by attaching the Strap Hook Fastener (Figure 1, number 6) directly to the Exterior Pad Fabric (Figure 1, number 1) or onto the strap itself. See Figure 3A.
3. To remove wrap, detach and unwrap Elastic Fastening Straps.

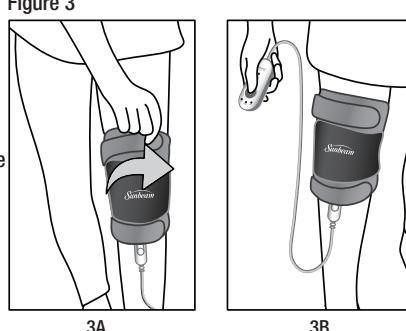
### HOW TO USE YOUR JOINT WRAP FOR HEAT THERAPY

#### USE ELECTRIC HEAT THERAPY FOR CONSISTENT, LONG LASTING, HIGH-LEVEL HEAT THERAPY

1. Ensure that the Gel Pack (Figure 1, number 9) is not inserted in the Mesh Gel Pack Pouch (Figure 1, number 3).
2. Make sure Cord Connector (Figure 4, letter A) is securely attached to the Pad Connector (Figure 4, letter B) by pushing both connector parts together until the Connector Release Tab (Figure 4, letter C) clicks into place for a tight fit.
- NOTE: Controller must be connected to heating pad before plugging unit into outlet.
3. After controller is connected to pad, plug into 120V AC electrical outlet.
4. To turn the pad on, press the Control Button (Figure 5, letter B).
5. Select from 3 heat settings: low, medium, or high, by repeatedly pressing the control button until the desired Heat Setting Indicator (Figure 5, letter C) lights up.
6. To turn the pad off, press the Control Button repeatedly until all of the Heat Setting Indicator lights turn off.


**CAUTION:**  
• DO NOT SIT ON, LEAN AGAINST, CRUSH OR FOLD PAD DURING USE.  
• DO NOT USE WHILE SLEEPING.  
• SEE ADDITIONAL WARNINGS IN THE IMPORTANT SAFETY INSTRUCTIONS SECTION.

Figure 4

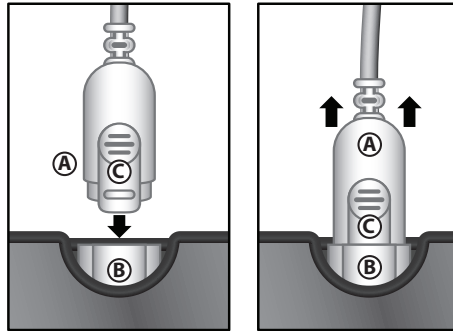A. Cord Connector  
B. Pad Connector  
C. Connector Release Tab

Figure 5

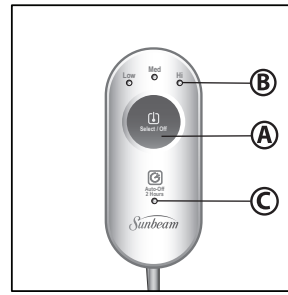A. Control Button  
B. Heat Setting Indicator  
C. Auto-Off Indicator

### AUTO-OFF FEATURE

Your heating pad controller includes an Auto-Off feature. This feature is designed for convenience and to conserve energy by automatically turning off after 2 hours of continuous use. The Auto-Off Indicator (Figure 5, letter C) will flash when the auto-off feature is activated and the pad shuts off. To restart, press the Control Button until the desired heat setting is reached.  
NOTE: After a power outage, the Auto-Off feature automatically resets the timer.

### FOR ON-THE-GO HEAT THERAPY USING GEL PACK

**BE EXTREMELY CAREFUL WHEN HANDLING PACK FROM MICROWAVE TO AVOID BURNS. HEAT ONLY IN A MICROWAVE. HEATING TIME MAY VARY DEPENDING ON MICROWAVE RATED OUTPUT.**

### Heating Instructions:

#### From Room Temperature:

1. Heat in microwave on HIGH in 30 second intervals until a safe and desired heat level is reached, with total heating time not to exceed 90 seconds. Massage and turn gel pack over in microwave between intervals to avoid hot spots.
2. Monitor gel pack in microwave while heating. If any swelling of the gel pack occurs, stop heating immediately and wait for the swelling to subside before removing.
3. Before removing gel pack from the microwave, carefully check for any signs of possible rupture or leakage. If any leakage is found, do not use the gel pack.

#### From Frozen:

1. Heat in microwave on HIGH for 60 seconds, then in intervals of 30 seconds until a safe and desired heat level is reached, with total heating time not to exceed 120 seconds. Massage and turn gel pack over in microwave between intervals to avoid hot spots.
2. Monitor gel pack in microwave while heating. If any swelling of the gel pack occurs, stop heating immediately and wait for the swelling to subside before removing.
3. Before removing gel pack from the microwave, carefully check for any signs of possible rupture or leakage. If any leakage is found, do not use the gel pack.

**INSERT GEL PACK IN PAD'S MESH POUCH PRIOR TO USE. NEVER APPLY DIRECTLY TO SKIN.**

#### CAUTION:

1. Consult a physician prior to use.
2. People with poor circulation, diabetes, nerve damage, paralysis or sensitive skin should not use this product.
3. Although non-toxic, discard gel pack if product leaks or is punctured.
4. Adult supervision is required for use on invalids or children. Do not use on infants or children under the age of 2 years.
5. Do not use the gel pack in conjunction with liniments, sports creams, balms or lotions.
6. Do not use on broken skin.
7. If gel pack is too hot, allow to cool before use.
8. Do not fall asleep with the gel pack.

Figure 6

### HOW TO USE YOUR JOINT WRAP FOR COLD THERAPY – SEE FIGURE 6

1. Place Gel Pack (Figure 1, number 9) in the freezer for at least 1 hour prior to use.
2. Disconnect controller from the pad by pushing down on the Connector Release Tab (Figure 4, letter C) and pulling the Cord Connector (Figure 4, letter A) away from the Pad Connector (Figure 4, letter B).  
Note: Heat therapy should never be applied while the gel pack is in use for cold therapy.
3. Once gel pack is cold, insert into the Mesh Gel Pack Pouch (Figure 1, number 3) and secure in place with the Gel Pack Pouch Fastener (Figure 1, number 4).  
Note: Always use product with the provided cover. Never apply directly to skin.
4. Apply wrap to the elbow or knee using the application instructions.

### STORAGE INSTRUCTIONS – SEE FIGURE 7

When FlexTemp™ Hot + Cold Joint Wrap is not in use, easily store away until next treatment.

1. When not in use, the Gel Pack can be stored in the Mesh Gel Pack Pouch.  
Note: For increased shelf life, store Gel Pack in cool, dry, and ventilated area away from incompatible substances. Keep away from sources of ignition. Do not store in freezer.
2. Neatly wrap cord and controller and place inside the Fabric Sleeve. Tightly wrap Elastic Fastening Straps around the pad and secure for compact, clean storage.

### WASHING INSTRUCTIONS

**CAUTION: DO NOT DRY CLEAN! DO NOT BLEACH! DO NOT USE WRINGER!**

1. Disconnect pad from outlet. DO NOT PULL BY THE CORD.
2. Disconnect Controller Cord from pad by pressing the tab on the Cord Side Connector and pulling on the Pad Side Connector.
3. Fold Elastic Fastening Straps and attach Strap Hook Fasteners to the Exterior Pad Fabric before placing in the washing machine.
4. Machine wash pad in cold water on gentle cycle or hand wash. Tumble dry on low heat and remove promptly from dryer or drape to dry.
5. DO NOT plug in pad and pad connector until completely dry.
6. DO NOT wash with the Gel Pack inserted.
7. DO NOT iron your heating pad.
8. For information on service, visit our website at [www.sunbeam.com](http://www.sunbeam.com).

### TROUBLE SHOOTING

| Problem                                                                                                             | Solution                                                                                                                                                                                                                                                                                                                                                                                                                                                                                    |
|---------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <ul style="list-style-type: none"><li>• No Heat</li><li>• No Power</li></ul>                                        | <ol style="list-style-type: none"><li>1. Verify tight connection between pad and controller.</li><li>2. Verify AC outlet is functional.</li><li>3. Verify no physical damage to the cord, controller or pad.</li></ol> <p>If failure continues, discontinue use. See warranty information to contact Sunbeam Consumer Service.</p>                                                                                                                                                          |
| <ul style="list-style-type: none"><li>• Pad has No Heat &amp; Controller 'Auto-Off' Indicator is Flashing</li></ul> | <p>Auto-Off feature is set to turn the unit off after 2 hours of continuous use.</p> <ul style="list-style-type: none"><li>• To re-start unit, press the Control Button until desired heat setting is reached.</li></ul>                                                                                                                                                                                                                                                                    |
| <ul style="list-style-type: none"><li>• Pad has a Hot Spot</li></ul>                                                | <ul style="list-style-type: none"><li>• Ensure the pad is not folded, bunched, or crushed in any area.</li></ul>                                                                                                                                                                                                                                                                                                                                                                            |
| <ul style="list-style-type: none"><li>• 'High' Heat Setting Indicator is Flashing</li></ul>                         | <p>The 'High' Heat Setting Indicator light will start flashing when connection is broken between pad &amp; controller while using the product.</p> <ol style="list-style-type: none"><li>1. To reset, simply connect pad &amp; controller together and press Control Button until desired heat setting is reached.</li><li>2. If flashing light persists, an error has occurred. Discontinue use of the product and see warranty information to contact Sunbeam Customer Service.</li></ol> |
| <ul style="list-style-type: none"><li>• Gel Pack is Leaking or Punctured</li></ul>                                  | <ul style="list-style-type: none"><li>• Discontinue use and dispose of gel pack.</li></ul>                                                                                                                                                                                                                                                                                                                                                                                                  |

### FCC CERTIFIED

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Note:** This product has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This product generates, uses and can radiate radio frequency energy and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the product and receiver.
- Connect the product into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

### 2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc. ("Sunbeam") warrants that for a period of two (2) years from the date of initial retail purchase, this product will be free from defects in material and workmanship, when used in accordance with the instructions provided. Sunbeam warrants the gel pack supplied with the unit will be free from defects in material and workmanship for one (1) year from the date of initial retail purchase.

Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period.  
Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty. Do NOT attempt to repair or adjust any electrical or mechanical functions on this product. Doing so will void this warranty. This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or any way change the terms and conditions of this warranty.  
This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use on improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

#### What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any express, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty. Sunbeam disclaims all other warranties, conditions or representations, express, implied, statutory or otherwise. Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party. Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you. This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

#### How to Obtain Warranty Service

In the U.S.A.  
If you have any question regarding this warranty or would like to obtain warranty service, please call 1 800 435-1250 and a convenient service center address will be provided to you.  
In the U.S.A., this warranty is offered by Sunbeam Products, Inc. located in Boca Raton, Florida 33431.  
**PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.**

For inquiries regarding recycling and proper disposal of this product, please contact your local waste management facility.

For product questions:  
Sunbeam Consumer Service  
USA: 1.800.435.1250  
[www.sunbeam.com](http://www.sunbeam.com)

© 2019 Sunbeam Products, Inc. All rights reserved. Distributed by Sunbeam Products, Inc., Boca Raton, Florida 33431. Sunbeam Products, Inc. is a subsidiary of Newell Brands Inc. (NASDAQ: NWL).

FlexTemp™ Joint Wrap\_901-825-002\_19ESM2 (US)  
Printed in China

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