
IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: To reduce risks of burns, electric shock, fire, and personal injury, this product must be used in accordance with the following instructions:

1. READ ALL INSTRUCTIONS
2. DO NOT USE WHILE SLEEPING.
3. NO UTILIZAR EN UN INFANTE OR ON AN ANIMAL.
4. THIS PAD IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, A PERSON WITH DIABETES, OR A PERSON WITH POOR BLOOD CIRCULATION.
5. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
6. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PAD FREQUENTLY.
7. NEVER USE PAD WITHOUT COVER IN PLACE.
8. **DO NOT USE IN OXYGEN ENRICHED ENVIRONMENT**
9. **DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.**
10. **DO NOT SIT ON OR AGAINST OR CRUSH PAD; AVOID SHARP FOLDS.**
11. **NEVER PULL THIS PAD BY THE SUPPLY CORD.**
12. **DO NOT USE THE CORD AS A HANDLE.**
13. CAREFULLY EXAMINE INNER COVER BEFORE EACH USE. DISCARD THE PAD IF INNER COVERING SHOWS ANY SIGNS OF DETERIORATION.
14. THIS UNIT SHOULD NOT BE USED BY CHILDREN WITHOUT ADULT SUPERVISION.
15. USE THIS PAD ONLY ON A 110-120 VOLT AC CIRCUIT. UNPLUG WHEN NOT IN USE.
16. DO NOT TAMPER WITH THIS PAD IN ANY WAY.
17. THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PAD DOES NOT FUNCTION SATISFACTORILY, SEE WARRANTY FOR CONSUMER SERVICE INFORMATION.
18. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS, ESPECIALLY ONES THAT CONTAIN HEAT-PRODUCING INGREDIENTS. SKIN BURNS COULD RESULT.
19. PLACE PAD ON TOP OF AND NOT UNDER THE PART OF BODY NEEDING HEAT.
20. SAVE THESE INSTRUCTIONS.

WARNING: Do not massage abdomen when pain is present. Do not massage areas that are swollen, inflamed, or covered with skin eruptions. Do not use on unexplained calf pain. Consult a physician before using any equipment of this type.

INCLUDES A UL/CUL LISTED TRANSFORMER
ELECTRIC HEATING PAD – FOR HOUSEHOLD USE ONLY.

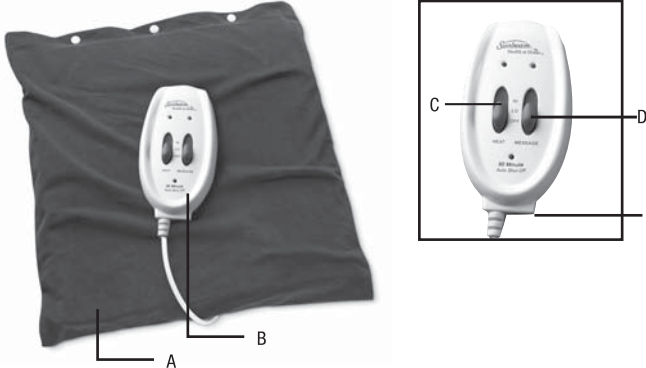Features

Your Sunbeam® Heat Plus Massage Heating Pad comes with the following:

- Soft cover with snap closure that is removable and machine washable.
- Slide control with 2 heat and 2 massage settings.

Features Of Your Heat Plus Massage Heating Pad

- Soothing warmth and vibration for sore muscles and temporary arthritic pain
- Use heat and massage therapy independently or together
- 2 heat and 2 massage levels
- Slide control with ON indicator light
- Auto shut-off after one hour

How to Use the Heat Plus Massage Heating Pad

Your Sunbeam® Heat Plus Massage heating pad combines soothing heat with vibration for sore muscles and temporary arthritic pain. The heat and massage features can be used together or independently.

For Dry Heat Application

- Lay the cloth cover on a flat surface with the embroidered side up. Open the three snaps.
- Insert the pad, with the heated side up, into the cover completely, fitting one side at a time.
- Make sure the cord and controller are on left side, and both remain on the outside of the cover.
- The pad will lie flat, filling the entire cover. Close the double-sided snaps on the pad to the matching snaps on the cover. If the snaps are not aligned, flip over the pad.
- Plug into nearest power supply.
- Plug power adapter into controller at plug port (F).
- Turn the controller on to desired heat and/or massage setting. Allow 15 minutes for the pad to reach the high heat setting.
- Relax and enjoy. **Do not use for more than 60 minutes at a time** (not counting warm-up time).
- When finished, turn controller off and unplug the heating pad.

CAUTION: Only use heating pad when inserted into cloth cover.

Where to Use the Heat Plus Massage Heating Pad:

Apply the heating pad to the area in need of treatment. Only use heating pad with cover in place.

Fabric Cover Washing Instructions:

Remove cover. Machine or hand wash; use a mild detergent. Wash with cool water only. Hang dry or machine dry on gentle, cool cycle.

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Trouble-Shooting:

Problem	Solutions
• Control Blinks • No Heat • No Power	• Unplug unit from wall outlet. • Verify no physical damage to the cord, controller or pad. • Verify tight connection between pad and controller. • Verify AC outlet is functional. • If blinking or failure continues, discontinue use. See Warranty Information to contact Consumer Services.
• Control Shuts Off	• Auto Off feature is set to turn unit off after one hour of use. • Re-start unit by sliding control button. • Verify tight connection between pad and controller.
• Pad has a Hot Spot	• Make sure pad is not folded or bunched up in any area
• Pad doesn't get hot enough	• Re-set Heat Selection by sliding control button. • Verify tight connection between pad and controller. • Heating Pad is designed to meet maximum heat output according to UL Standards.

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2-YEAR LIMITED WARRANTY

Sunbeam Products, Inc., or if in Canada, Sunbeam Corporation (Canada) Limited (collectively "Sunbeam") warrants that for a period of two years from the date of purchase, this product will be free from defects in material and workmanship. Sunbeam, at its option, will repair or replace this product or any component of the product found to be defective during the warranty period. Replacement will be made with a new or remanufactured product or component. If the product is no longer available, replacement may be made with a similar product of equal or greater value. This is your exclusive warranty.

This warranty is valid for the original retail purchaser from the date of initial retail purchase and is not transferable. Keep the original sales receipt. Proof of purchase is required to obtain warranty performance. Sunbeam dealers, service centers, or retail stores selling Sunbeam products do not have the right to alter, modify or any way change the terms and conditions of this warranty.

This warranty does not cover normal wear of parts or damage resulting from any of the following: negligent use or misuse of the product, use of improper voltage or current, use contrary to the operating instructions, disassembly, repair or alteration by anyone other than Sunbeam or an authorized Sunbeam service center. Further, the warranty does not cover: Acts of God, such as fire, flood, hurricanes and tornadoes.

What are the limits on Sunbeam's Liability?

Sunbeam shall not be liable for any incidental or consequential damages caused by the breach of any expressed, implied or statutory warranty or condition. Except to the extent prohibited by applicable law, any implied warranty or condition of merchantability or fitness for a particular purpose is limited in duration to the duration of the above warranty.

Sunbeam disclaims all other warranties, conditions or representations, expressed, implied, statutory or otherwise.

Sunbeam shall not be liable for any damages of any kind resulting from the purchase, use or misuse of, or inability to use the product including incidental, special, consequential or similar damages or loss of profits, or for any breach of contract, fundamental or otherwise, or for any claim brought against purchaser by any other party.

Some provinces, states or jurisdictions do not allow the exclusion or limitation of incidental or consequential damages or limitations on how long an implied warranty lasts, so the above limitations or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights that vary from province to province, state to state or jurisdiction to jurisdiction.

How to Obtain Warranty Service

If you have any question regarding this warranty or would like to obtain warranty service, please call a and a convenient service center address will be provided to you, or visit our web site at www.sunbeamhealth.com.

In the U.S.A. - This warranty is offered by Sunbeam Products, Inc., located in Boca Raton, Florida 33431, 1-800-435-1250.

In Canada, this warranty is offered by Sunbeam Corporation (Canada) Limited, located at 5975 Falbourn Street, Mississauga, Ontario L5R 3V8, 1-800-667-8623.

PLEASE DO NOT RETURN THIS PRODUCT TO ANY OF THESE ADDRESSES OR TO THE PLACE OF PURCHASE.

Printed in China

© 2005 Sunbeam Products, Inc. doing business as: Jarden Consumer Solutions. All rights reserved. In Canada: Imported and distributed by Sunbeam Corporation (Canada) Limited doing business as: Jarden Consumer Solutions, 5975 Falbourn Street, Mississauga, Ontario L5R 3V8.

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