

## Medical Materials & Technologies

# 3M™ Microfluidic Diagnostic Tape, 9965 Double-Sided White Tape




### Product Description:

3M™ Microfluidic Diagnostic Tape 9965, consists of a white polyester film coated on both sides with a pressure sensitive acrylate adhesive. The tape is supplied between two clear, silicone coated polyester release liners. The product is designed for use in diagnostic medical devices such as test strips and lab-on-a-chip.

### Features and Benefits:

- Hydrophobic adhesive resists clouding
- Adheres to a variety of substrates
- White polyester carrier
- Double sided adhesive coating
- Non-hemolytic
- Non-toxic to mammalian cells
- Low extractables
- Designed for lab-on-a-chip applications

### Composition:

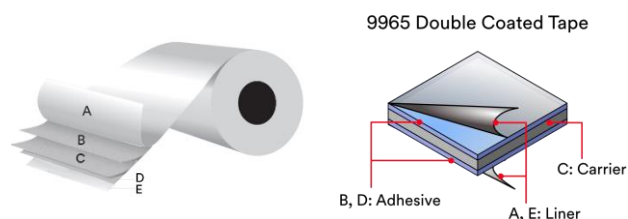

|                  |   |
|------------------|---|
| A: Release Liner | 2 mil (51 um) Clear silicone coated polyester release liner |
| B: Adhesive      | 0.7 mil (18 um) Acrylate adhesive                           |
| C: Carrier       | 2 mil (51 um) White polyester film                          |
| D: Adhesive      | 0.7 mil (18 um) Acrylate adhesive                           |
| E: Release Liner | 2 mil (51 um) Clear silicone coated polyester release liner |

### Roll Description:

|                  |  |
|------------------|--|
| Roll Length      | Maximum 350 yd (320 m) ± 2%                |
| Roll Width       | Maximum 48 in (122 cm) ± 0.06 in (1.52 mm) |
| Splices per Roll | Maximum 3 per 350 yd (320 m)               |

Note: Default 6 in (152 mm) plastic core

### Adhesion & Liner Release:

|  |                       |
|--|-----------------------|
| Adhesion to steel, 180° Peel*          | 21.1 oz/in (2.3 N/cm) |
| Adhesion to polyester, 180° Peel*      | 34.1 oz/in (3.7 N/cm) |
| Adhesion to polycarbonate, 180° Peel*  | 34.9 oz/in (3.8 N/cm) |
| Adhesion to acrylic, 180° Peel*        | 38.1 oz/in (4.2 N/cm) |
| Adhesion to polypropylene, 180° Peel*  | 18.9 oz/in (2.1 N/cm) |
| Liner Release moderate side, 180° Peel | 47.9 gm/in (0.2 N/cm) |
| Liner Release easy side, 180° Peel     | 8.6 gm/in (33 mN/cm)  |

Represents typical value. \*Adhesion tested after 5-minute dwell time.

**Biocompatibility Testing:**

The adhesive used in product 9965, in conjunction with a similar product with a different backing, has been subjected to the following safety evaluations.

| Study Type            | Test Method               | Standard              | GLP | Results                     |
|-----------------------|---------------------------|-----------------------|-----|-----------------------------|
| In vitro cytotoxicity | Cell lysis & cytotoxicity | ISO 10993-5: $\leq 2$ | Yes | Reactivity grade: 0         |
| MEM elution           | Cell lysis & cytotoxicity | ISO 10993-5: $\leq 2$ | Yes | Reactivity grade: 0         |
| Hemolytic Index       | Interaction with Blood    | ISO 10993-4: $< 2\%$  | Yes | Non-Hemolytic (index: 1.2%) |
|                       | Test article extract      | ISO 10993-4: $< 2\%$  | Yes | Non-Hemolytic (index: 0.0%) |

**Recommended Storage Conditions & Shelf Life:**

Product as supplied in original packaging will maintain certified properties for a period of two years from date of manufacture when stored at room temperature 68-77 °F (20-25 °C).

**Product and Safety Information:** User is solely responsible for determining the suitability of 3M samples and products for the intended use including any necessary safety or toxicity assessment. 3M will provide Material Safety Data Sheets or equivalent and summary results of biocompatibility testing upon request. In every case before using any product in full scale production users should conduct their own tests to determine to their own satisfaction whether the product is of acceptable quality and is suitable for their particular purposes under their own operating conditions.

**Notice:** Nothing contained herein shall be construed to imply the nonexistence of any relevant patents or to constitute a permission, inducement or recommendations to practice any invention covered by any patent, without authority from the owners of this patent.

#### Warranty Information

All statements, technical information and recommendations herein are based on tests 3M believes to be reliable, but the accuracy or completeness thereof is not guaranteed. 3M warrants only that products will meet 3M's specifications at the time of shipment to the customer. 3M does not offer any other warranty and does not warrant the performance, safety, or such other characteristics of products in combination with other materials. 3M specifically DOES NOT warrant products for any intended or unintended uses (whether or not foreseeable); for compatibility or suitability with other components or compatibility with any methods of manufacture or conversion. The foregoing warranty is made in lieu of all other warranties, expressed or implied, including the implied warranties of merchantability, fitness for a particular purpose and freedom from non-infringement.

**Limitation of Remedies:** If products are proven not to meet 3M's specifications, the sole and exclusive remedy available and 3M's only obligation shall be, at 3M's option, to replace such quantity of Products which are proven out of specification or to refund the purchase price paid for Products.

**Limitations of Liabilities:** The remedies provided herein are exclusive remedies against 3M for any alleged or actual nonconformance to specifications or defect or other failure in products or for 3M's performance of its supply obligations. Under no circumstances is 3M liable for any direct, indirect, incidental, special, or consequential damages (including lost profits) in any way related to the product under any theory of law including, but not limited to, negligence and strict liability.

#### Ordering Information

To place an order, call our Customer Service at: **800-742-1994 (U.S.)**. For all orders outside the United States, please contact your local country representative or check out **3M.com/MedTechFindaRep** for your region's contact information.

To have a sales representative contact you, to request samples or clinical and safety summaries, please contact us at 3M HELPLINE 800-228-3957 (U.S.) or check out **3M.com/MedTechContactUs** for your region's contact information.

Visit our website: **www.3M.com/MedTech** for product and services information, news, events, new product highlights or to make a direct inquiry.



#### Medical Materials & Technologies

3M Center, Building 275-5W-05  
St. Paul, MN 55144-1000 USA

Phone 1-800-228-3957  
Web 3M.com/MedTech

3M is a trademark of 3M.

Please recycle. Printed in U.S.A.

© 3M 2022 All rights reserved.

Effective Sep 2022