
Opticlude™

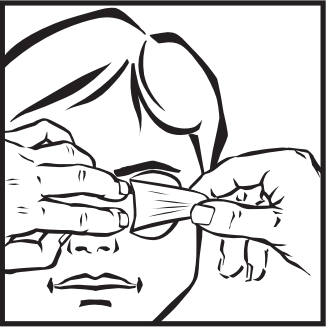Removeal

- At the nose side of eye patch, grasp and loosen the outer edge of the eye patch.

- Support the skin with your fingers.

Remove the tape low and slow, peeling the eye patch towards the ear, keeping it close to the skin surface and pulled back over itself. Please be aware that removing the eye patch at an angle will pull the skin surface and increase the risk of redness and irritation.

- As the patch is removed, continue to support the newly exposed skin.

MD

CE

Made in Poland

3M Deutschland GmbH
Health Care Business
Carli-Schurz-Str.1
41453 Neuss, Germany

©2015 3M. All rights reserved. 3M and Opticlude are trademarks of 3M. 3M.com DV-01063-0163-7/0819 3100021585/0819

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Opticlude™

Orthoptic Eye Patch for occlusion therapy

Description: 3M™ Opticlude™ Orthoptic Eye Patch consists of a soft nonadherent pad bonded to a skin-friendly breathable and conformable backing coated with a hypoallergenic adhesive. The dressing is breathable and conforms to body contours.

Intended Use: 3M™ Opticlude™ Orthoptic Eye Patch is used in eye occlusion treatment of Amblyopia.

Precautions: Follow directions provided by your eye care professional regarding use of 3M™ Opticlude™ Orthoptic Eye Patch. Use only on intact skin. Follow instruction for use to minimize risk of skin redness and irritation. In case of any skin reactions, immediately stop using the patch and consult your doctor. For patients with sensitive skin we recommend to use Opticlude™ Silicone eyepatches. This is a single use device. Reuse may result in compromising product integrity or lead to device failure.

Instructions for use:

- Application**
- For best adhesion, first clean and dry skin. Remove paper liner from the eye patch.
 - Gently close both eyes keeping facial muscles relaxed. Wrinkles around the eye should be avoided.
 - Apply 3M™ Opticlude™ Orthoptic Eye Patch with the narrow end towards the nose. Do not stretch the patch or pull skin.
 - Secure the pressure sensitive adhesive tape by applying light pressure with fingertips around the outside of the eye patch.

- Removeal**
- At the nose side of eye patch, grasp and loosen the outer edge of the eye patch.
 - Support the skin with your fingers. Remove the tape low and slow, peeling the eye patch towards the ear, keeping it close to the skin surface and pulled back over itself. Please be aware that removing the eye patch at an angle will pull the skin surface and increase the risk of redness and irritation.
 - As the patch is removed, continue to support the newly exposed skin.

Storage/ Shelf Life/Disposal
The device can be stored under normal conditions. The device can be disposed as normal household waste. In case it is contaminated, please follow your local waste disposal regulations. After the shelf life has expired, the device shall not be used but disposed as indicated.

How supplied / Ordering Information
Opticlude Orthoptic Eye Patches are available in multiple sizes and packskises with a variety of different designs. Please contact your local 3M subsidiary for details about offerings in your country.

Note: Please report a serious incident occurring in relation to the device to 3M and the local competent authority (EU) or local regulatory authority.

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