

Panda™ Warmer

Care that can help you deliver
seamless care for your tiniest patients



GE HealthCare

Let's give every baby the strongest start we can

You want to give every baby, no matter how small the size or big the challenge, the chance for the best start in life. But when babies need a little extra care and attention and can't be with mom, the most important tools you have are your clinical know-how paired with your empathy to give newborns the chance they deserve.

Panda Warmers are an all-in-one integrated solution that you can use to help babies go home healthy. We designed them to create a more inclusive and collaborative environment for care providers and their patients. In addition, Panda Warmers provide integrated technology solutions to help you deliver responsive, appropriate care.



Integration that allows you to focus on your most important task

As a labor and delivery provider, you're a master at juggling multiple tasks while staying vigilant for any number of possible challenges. Panda Warmers put all your tools in one place, helping you focus on your most important task—caring for newborn patients and their moms.

Simple



A bright idea

The standard observation light lets you adjust light intensity so you can work without creating unnecessary discomfort for the newborn. The optional integrated procedure light lets you easily focus the light where needed.



Brilliant, full-color controls

The control panel & full-color display are accessible and easy to use. Patient information is displayed vividly on the large screen.



Easily expanded

Our vertical dovetail dual rail system is incorporated for easy addition of accessories and ancillary equipment to simplify your equipment management.

Efficient



Warm-up with no nuisance alarms

Automatic warm-up mode helps ensure the bed is warm when the newborn arrives. By eliminating nuisance alarms, we help keep the environment quiet and stress-free for caregivers, babies, and parents.



Heat where you need it

Our innovative Hourglass Heat Profile bathes the baby in soothing warmth with uniform coverage across the mattress, while helping keep the caregiver cool and comfortable.

Intelligent



Hands-free alarm silence

Maintain a clean environment and keep the area quiet and soothing by simply moving your hand in front of the sensor to silence the alarm.



Accurate temperature control

Powered by the same closed-bed, heat-controlled algorithm as our trusted Giraffe™ OmniBed™ Carestation™.


Easily customize your Panda Warmer

Make the most of your investment and provide the best possible care with these optional features.



Integrated resuscitation

Available with minimal setup to meet the Neonatal Resuscitation Program (NRP) and International Liaison Committee on Resuscitation (ILCOR) guidelines, this feature helps standardize resuscitation protocols across the perinatal care area. Integrated gas cylinder holders and easily-viewable tank pressure gauges enable smooth workflow and transport convenience. Integrated t-piece system with a built-in manometer and integrated air/oxygen blender keeps important tools at your fingertips.

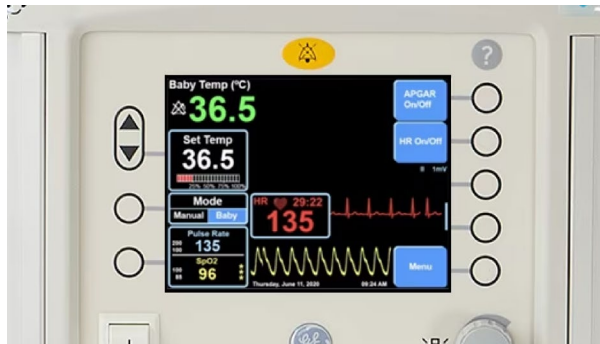


Integrated SpO₂ and ECG heart rate

The integrated ResusView™ Solution*, available with Panda iRes warmers, provides accurate ECG heart rate and SpO₂ blood oxygen levels within seconds, allowing you to provide timely and appropriate interventions.



Uninterruptible Power Supply (UPS)*

The UPS battery keeps the baby safe and secure from power interruptions for 15 minutes.**

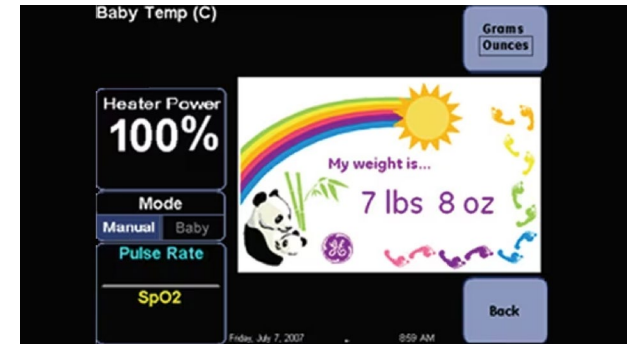

Integrated scale

Baby's first weight is available right at the bedside. No need to move the baby to a cold scale.



X-ray tray

With an integrated X-ray cassette tray, baby can undergo a radiographic exam right in the warmer with uninterrupted thermal support.

Clinical confidence right after birth

You give the world's newborns a strong start at life, and we're here to support you with reliable, advanced technology that's also incredibly easy to use. Our Panda Warmer ResusView solution, available with the Panda iRes Warmer or as an upgrade to other Panda Warmers, helps give your staff clinical confidence in those vital moments.

Heart rate in 6 seconds***

An accurate heart rate assessment immediately after birth is critical for guiding resuscitation efforts. With ResusView, ECG heart rate and blood oxygen saturation level (SpO2) are integrated into your warmer. As an added benefit, ResusView gives you a convenient view of this data with a continuous ECG heart rate number and waveform displayed on the Panda Warmer screen, something palpation or auscultation won't provide.





GE HealthCare

About GE HealthCare Technologies Inc.

GE HealthCare is a leading global medical technology, pharmaceutical diagnostics, and digital solutions innovator, dedicated to providing integrated solutions, services, and data analytics to make hospitals more efficient, clinicians more effective, therapies more precise, and patients healthier and happier. Serving patients and providers for more than 100 years, GE HealthCare is advancing personalized, connected, and compassionate care, while simplifying the patient's journey across the care pathway. Together our Imaging, Ultrasound, Patient Care Solutions, and Pharmaceutical Diagnostics businesses help improve patient care from diagnosis, to therapy, to monitoring. We are a \$19.6 billion business with 51,000 colleagues working to create a world where healthcare has no limits.

Follow us on [LinkedIn](#), [X](#) (formerly Twitter), and [Insights](#) for the latest news, or visit our website <https://www.gehealthcare.com/> for more information.

Disclaimers

* Not available in all markets

** Dependent on load conditions

*** 60-119 bpm

GE HealthCare reserves the right to make changes in specifications and features shown herein, or discontinue the product described at any time without notice or obligation. Contact your GE HealthCare representative for the most current information.

© 2025 GE HealthCare. Panda, Giraffe, OmniBed, Carestation, and ResusView are trademarks of GE HealthCare.
GE is a trademark of General Electric Company used under trademark license.

August 2025
JB29324XX_2