

The bulk fill flowable that's easier to inject

3M™ Filtek™ Bulk Fill Flowable Restorative is available in a new ergonomic syringe designed for improved comfort and easier injecting. This gives you more control and visibility for precise placement in deep restorations.

Reduce bubbles and waste

Virtually no bubbles or material 'run-on' during dispensing gives you better control



Bendable cannula


Easy-to-bend without kinking for better access

No more guessing

The green syringe barrel indicates the remaining material volume. When your syringe is empty, only the white plunger shows.

Easy to hold and inject

Triangle finger plate and plunger make it quick and easy to personalise your grip

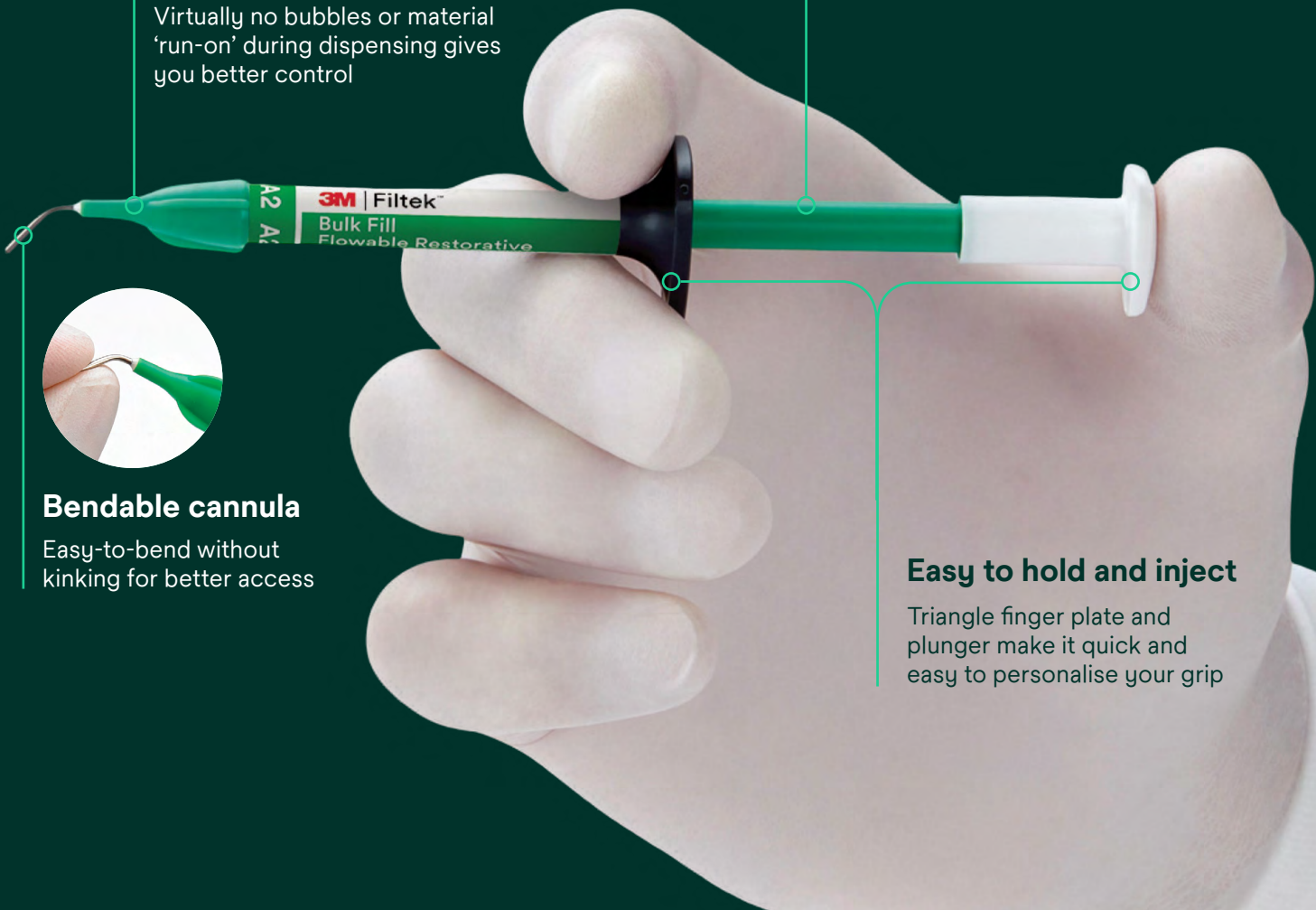

Key benefits

- Virtually no bubbles or material 'run-on'
- Easy-to-bend cannula for improved access
- Ergonomic design is easy to hold and inject
- Injectable, one-step placement up to 4 mm
- Excellent adaptation and low shrinkage stress

Can be warmed for injection molding[^]

Warm up to 70°C (158°F) for up to one hour, up to 25 times.



Posterior restoration made simple



Prepare the tooth and apply an adhesive e.g. 3M™ Scotchbond™ Universal Adhesive.



Dispense up to 4 mm of 3M™ Filtek™ Bulk Fill Flowable Restorative into the prepared cavity and cure, e.g. with 3M™ Elipar™ DeepCure-S LED Curing Light.



Place a 2 mm (minimum) capping layer of an aesthetic composite, such as a 3M™ Filtek™ Universal Restorative, and cure.

Clinical photos courtesy of Dr. Reich, Seefeld, Germany.

[^] Data on file.

1. Including small, non stress-bearing occlusal restorations.
2. Repair of acrylic temporary and composite resin.
3. Where at least half the coronal tooth structure is remaining to provide structural support for the crown.

Indications


(Incremental up to 4 mm)



Base
(Class I – II)



Liner
(Class I – II)

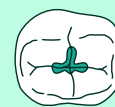

Minimally Invasive
Dentistry (MID)¹



Class III



Class V



Repair of indirect
(small defects)



Repair²



Undercut
block-out



Pit & fissure
sealant



Core Build-up³

Refer to Instructions for Use for full list of indications.

Ordering information



3M™ Filtek™ Bulk Fill Flowable Restorative Refills

Syringe refills include:
2 – 2 g syringes; 20 dispensing tips.

Shade	Item no.
A1	4862A1
A2	4862A2
A3	4862A3
Universal	4862U
Item no.	Product info.
3740T	3M™ Filtek™ Bulk Fill Flowable Restorative Dispensing Tips, 19-gauge, Green (20)



Solventum Dental Solutions

KCI Medical Australia Pty Ltd
Level 3, Building A, 1 Rivett North Ryde, NSW 2113
1300 363 454

KCI New Zealand Unlimited
Suite 1701, Level 17, PwC Tower, 15 Customs Street West
Auckland Central, Auckland 1010
0800 808 182

www.solventum.com

Please always refer to the product Instructions for Use.

© Solventum 2026. Solventum, the S logo, Elipar and Filtek are trademarks of Solventum or its affiliates. 3M, the 3M logo and Scotchbond are trademarks of 3M. All other trademarks are owned by their respective owners.