

## 3M™ Tube Securement Device

This strong, easy-to-use solution secures non-critical medical tubes made from silicone, PVC and silicone coated latex, with two sizes spanning 6-24 French. It is gentle to skin and designed to support pressure injury and medical adhesive-related skin injury (MARSII) prevention practices.


In addition to urinary catheters, the product is indicated for non-critical tube securement or in conjunction with sutures and/or primary securement dressings per facility protocol for surgical drains and vascular access applications.

### Uniquely shaped

Designed to lift tubes off skin and help promote consistent application.

### Sterile device

In single-use packaging.

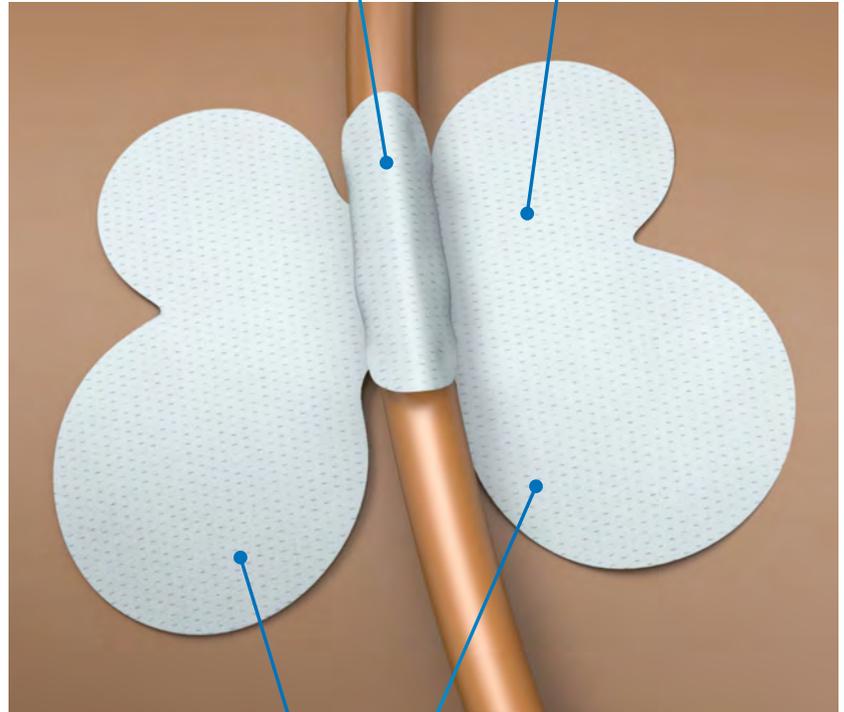

### Reinforced soft cloth backing

Provides extra strength and flexibility.

### Two sizes

Small device for medical tubes and catheters 6–12 Fr and large device for medical tubes and catheters 12–24 Fr.



**Skin performance:** The feather designates products that deliver the securement power you need while minimizing damage to skin.

### Reliable securement

Unique adhesive secures medical tubes and catheters for up to 4 days.<sup>1</sup>

### Ease of use

Easy to apply and removes cleanly in one piece; allows for standardization of technique.

### Skin safe

Non-sensitizing silicone adhesive has low potential to cause dermatitis with minimal skin damage and pain upon removal.<sup>1</sup>

### Designed to help reduce pressure injury risk

With no hard components to cause pressure, product's unique shape minimizes contact between skin and tube.


When it comes to tube securement selection there are a variety of factors to keep in mind, including standardization, securement performance, skin protection and clinician ease of use. Consider the results of comparative studies, such as the one below, to help inform your decision.

### A comparison of the 3M™ Tube Securement Device to the StatLock® Foley Stabilization Device.<sup>1</sup>

#### Wear time:

Dependability and length of wear are important in choosing the right device.


The 3M™ Tube Securement Device had a statistically higher percentage of devices with the tube in place after a 72 hour wear than StatLock® Foley Stabilization Device.\*

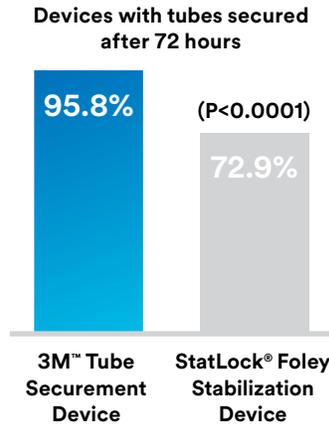

#### Skin integrity:

Skin is the body's first line of defense, so it is important to protect it.

Consider this comparison of the impact on the skin's barrier function (measured by average change in Transepidermal Water Loss (TEWL) from pre-treatment).

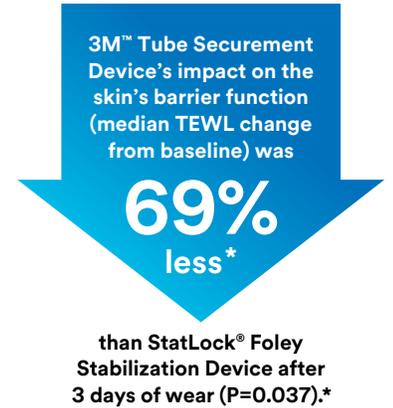

#### Removal time:

Efficient removal is critical when clinician time is limited.


The average removal time for 3M™ Tube Securement Device was 13 seconds compared to 126 seconds for StatLock® Foley Stabilization Device (P<.0001).



#### Wipes for removal:

Extra supplies can add time and cost to a procedure.

On average, the StatLock® Foley Stabilization Device required 12 alcohol wipes for removal, while the 3M™ Tube Securement Device does not require a removal agent.\*

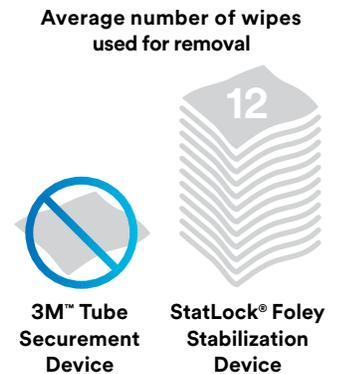

\*As tested in a comparative, prospective clinical study on 24 healthy human subjects after 3 days of wear followed by device removal from skin per IFU.

Feel secure in your choice. Trust the science behind 3M™ Tube Securement Device for reliable performance.

To order, contact your 3M Health Care Sales Representative or call the 3M Health Care Customer Helpline at 1-800-228-3957. Visit [3M.com/Medical](http://3M.com/Medical) for more information.

Refer to product Instructions for Use (IFU) for other important information.

#### Ordering Information

Product	Product Name	Product Number	Product Size	Tube/Catheter Size	Devices/Box	Boxes/Case
	3M™ Tube Securement Device, Large	1501U	2.2 in x 3.5 in (5,5 cm x 9,0 cm)	12-24 French	25	4
	3M™ Tube Securement Device, Small	1500U	1.7 in x 2.4 in (4,3 cm x 6,1 cm)	6-12 French	25	4



**3M Company**  
2510 Conway Avenue  
St. Paul, MN 55144-1000 USA

Phone 1-800-228-3957  
Web [3M.com/Medical](http://3M.com/Medical)

1. 3M Data on File. EM-05-014985.

© 2023 3M. All rights reserved.  
3M and the other marks shown are marks and/or registered marks. Unauthorized use prohibited.  
3M marks used under license in Canada.  
All other marks are property of their respective owners.  
GL\_ENG\_70-2013-1339-5